
No-Touch Forehead Thermometer Manual

Appearance & Functions

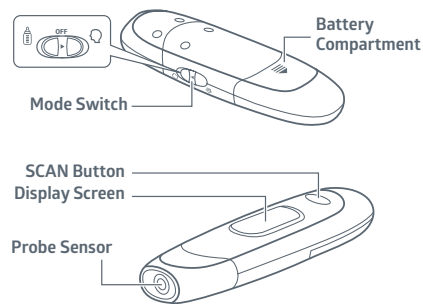Introduction


Thank you for choosing the infrared forehead thermometer, MDSNOTOUCH. This innovative medical device relies on advanced infrared (IR) technology to measure temperature instantly. The thermometer measures the infrared heat generated by the surface of the skin over the vessel and its surrounding tissue. The thermometer then converts it into a temperature value shown on the LCD screen.

Battery Installation


1. Remove battery cover.
2. Insert the included CR 2032 Battery with the + sign facing up. Reattach the battery cover.

LCD Screen Display Marks

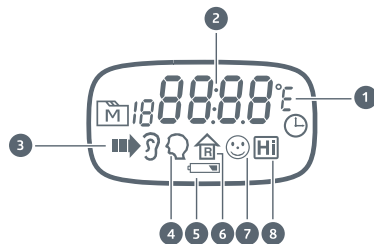
1. Temperature unit
2. Temperature display
3. Temperature scanning in progress
4. Forehead temperature indicator
5. Low battery indicator
6. Object surface temperature indicator
7. Smiling face indicator (Temperature falls within normal range)
8. High temperature indicator


Select Temperature

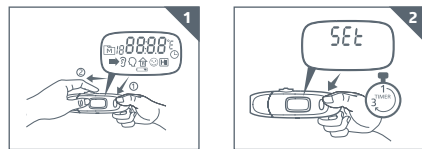1. Slide the Mode Switch from 'off' to 'forehead' mode and then immediately press and hold the 'Scan' button.
2. Keep holding down the 'Scan' button for 3-4 seconds until you see the word 'Set' appear on the screen.

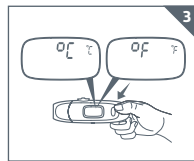
3. After the word 'Set' appears, press the 'Scan' button to toggle between C and F temperature units.

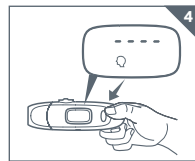4. After you have selected the preferred temperature unit, press and hold 'Scan' for 3-4 seconds until the screen readies for a forehead measurement (you will see the forehead symbol appear). The next measurement will display in the selected unit.


Tips for Performing A Temperature Measurement

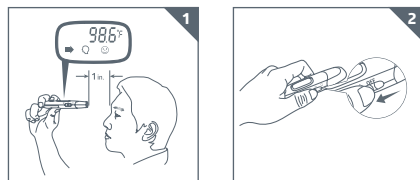1. To take a temperature, the thermometer must be held steady no more than 1 inch (2.3 cm) from the forehead. Sweat and other obstructions, such as hair, will affect the accuracy of the reading. Please make sure that the forehead is clear and dry before attempting a measurement. DO NOT hold the thermometer against the forehead.

2. When the measurement is complete, use the Mode Switch to turn the thermometer off.

Perform A Measurement

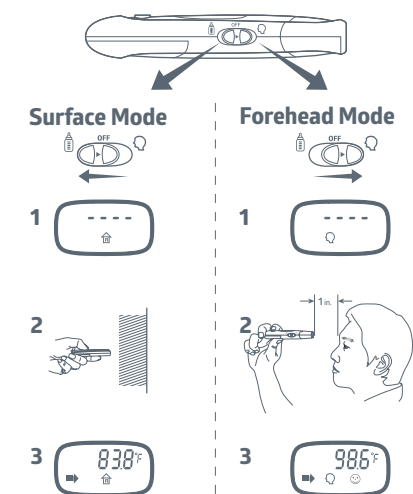
NOTE: While in surface mode, the measurement will perform continuously if you press and hold the scan button without release.

Result Indicators

	Temperature falls within measurement range: 93.9°F to 100.2°F (34.4°C to 37.9°C) and the LCD backlight will display "Green".
	Temperature falls within measurement range: 100.4°F to 108°F (38.0°C to 42.2°C) and the LCD backlight will display "Red".
	Temperature falls outside the displayed range: <ul style="list-style-type: none"> • Hi-temperature $\geq 108.1^\circ\text{F}$ (42.3°C). • Lo-temperature $\leq 93.7^\circ\text{F}$ (34.3°C) is shown on the display with "Blue" LCD backlight.

NOTE: Those indicators only appear for "Forehead Mode" measurements

Replacing The Battery

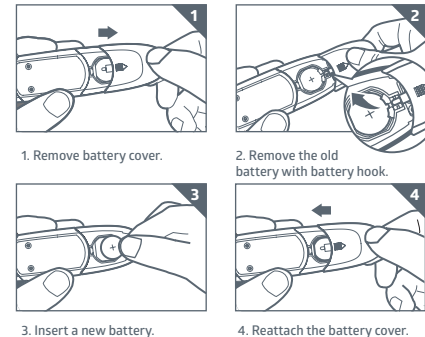Troubleshooting

If you get a "Lo" results indicator when measuring a surface temperature, you are in the wrong Mode.

Message	What it means	What to do
	Environmental temperature is below or above system operation range.	Put the thermometer under operating temperature range, 50°F to 104°F (10°C to 40°C)
	Problem with the thermometer.	Review the instructions and re-start the measurement. If the steps do not work, contact Medline customer service.
	Battery is low and "Lo" appears on LCD.	The meter is functional and the result remains accurate, but it is time to change the battery.
	The battery can't provide enough power for a test. After 3 beeps, the meter will turn off automatically.	Replace the battery immediately.
	Attempting to measure a surface temperature but are in "Forehead Mode".	Slide the Mode Switch to "Surface Mode" then repeat the measurement.

Detailed Temperature Information

Body temperature can vary from one individual/person to the next. It also varies by location on the body and time of day. Below shows the statistical normal ranges from different temperature sites. Please keep in mind that temperatures measured from different sites should not be directly compared, even if they are taken at the same time.

Table *1 Body Site Normal Temperature Range

Oral	1°F (0.6°C) or more above or below 98.6°F (37°C)
Rectal/Ear	0.5°F to 1°F (0.3°C to 0.6°C) higher than oral temperature
Axillary (armpit)	0.5°F to 1°F (0.3°C to 0.6°C) lower than oral temperature

*1. Body Temperature at WebMD; website: <http://firstaid.webmd.com/body-temperature> retrieved at 2010 Jan 7.

Warranty

Medline Industries, LP warrants this precision digital thermometer against any defects in material and workmanship for a period of one (1) year after the date of purchase. This guarantee excludes cases of freight damage, complete water immersion, tampering, clear abuse, misuse, or accidents. Medline Industries, LP will, at its discretion, repair or replace this thermometer during the warranty period without charge. Claims must be filed together with proof of purchase within this warranty period. Battery life and/or function is not included under this warranty. No representative or person is authorized to assume for us any liability in connection with the sale of the products of Medline Industries, LP. This warranty gives you specific legal rights, and you may also have other rights that may vary from state to state. Prior to shipping any product, please contact either Medline Industries, LP at 1-800-MEDLINE or an authorized Medline sales representative with any warranty concerns.

Important Safety Instructions


1. Close supervision is necessary when using this device around children or persons needing special care.
2. Use the thermometer only for the intended use described in this manual.
3. Do not use the thermometer if it is not working properly, or if it is damaged.
4. Do not use accessories that are not supplied or recommended by the manufacturer.
5. Proper maintenance is essential to the longevity of your device. If you are concerned about the accuracy of the device, please contact the manufacturer.

Cautions & Warnings

1. As with any thermometer, proper technique is crucial to getting accurate temperature readings. **Please read these instructions thoroughly and carefully before using the thermometer.**
2. Always operate the thermometer in an environment with a temperature range of 50°F to 104°F (10°C to 40°C), and relative humidity less than 95%.
3. Always store the thermometer in a cool and dry place: temperatures between -4°F to 140°F (-20°C to 60°C); relative humidity less than 95%.
4. **Avoid direct sunlight and avoid dropping the thermometer.**
5. **Do not** disassemble the thermometer.
6. Basic safety precautions should always be observed.
7. This thermometer is not intended to be a substitution for a consultation with your physician.

Cleaning & Disinfecting Procedures

1. Wipe clean the complete thermometer exterior surface including the display and buttons with a Micro-Kill+ disinfecting wipe (EPA Reg. No. 59894-10-37549). Do NOT use disinfecting wipe or tissue paper on the probe lens. Wipe the probe lens gently only with a dry cotton swab. The probe is NOT waterproof.

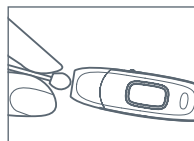2. Disinfect the meter in the same way as cleaning with another disinfecting wipe.
3. Keep the meter surface wet with disinfection solution for the contact time listed on the Micro-Kill+ disinfecting wipe's instructions for use (EPA Reg. No. 59894-10-37549).
4. Allow the thermometer surface to dry completely.
5. Discard the used wipes and never reuse them. Wash your hands thoroughly with soap and water after handling the thermometer to avoid contamination.

Improper system cleaning and disinfection may result in thermometer malfunction. If you notice any of the below signs after cleaning and disinfecting your thermometer, stop using the device and call Medline at 1-800-MEDLINE:

- Thermometer cannot be turned on
- LCD display cracks or becomes cloudy
- Buttons no longer function
- Thermometer outer casing cracks
- Color or paint/printing on housing is abnormal
- Scratches or abrasions on meter are higher than acceptable

This device has been validated to withstand 1500 cleaning and disinfection cycles using the recommended disinfecting wipe (EPA Reg. No. 59894-10-37549).

NOTE:

- We recommend that you clean the thermometer regularly to maintain accurate results.
- Do **NOT** clean and disinfect the thermometer while performing tests.
- Do **NOT** allow cleaning and disinfecting solution to get in the probe.
- If you do get moisture in the probe, wipe it away gently only with a dry cotton swab.
- Always dry the thermometer thoroughly before using it.
- Do **NOT** spray the thermometer directly with cleaning solutions especially those containing water (i.e. soapy water), as this could cause the solution to enter the case inside and damage the electronic components or circuitry.
- Do **NOT** use abrasive cleaners to clean the thermometer.

Symbol Information

-  Consult instructions for use
- IP22** Ingress protection rating
-  Manufacturer
- SN** Serial number
-  Caution, consult accompanying documents
-  Authorized representative in the European Community
-  Collection for electrical & electronic equipment
-  Type BF Equipment
-  Temperature limitation
-  Humidity limitation

Reference Standards

Device Standard:

Device Corresponds to the requirements of the standard for infrared thermometers

ASTM E1965-98, EN ISO 80601-2-56:2012,

EN 60601-1-2:2007/AC:2010,

EN 60601-1:2006+A12:2014, EN 60601-1-6:2010

Electromagnetic Compatibility:

Device fulfills the stipulations of the standard EN 60601-1-2

The stipulations of EU-Directive 93/42/EEC for Medical Device Class II devices have been fulfilled.

311-1115300-002

www.medline.com
Manufactured for: Medline Industries, LP, Three Lakes Drive, Northfield, IL 60093 USA. Made in Taiwan
1-800-MEDLINE V3 RD23TTH

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