

A photograph of a nurse in blue patterned scrubs sitting on a bed. A teal underpad is placed on the white bedsheet. The nurse's arms are crossed, and the underpad is positioned in front of her. The background is a plain, light-colored wall.

EFFECTIVE MOISTURE MANAGEMENT.

PerforMAX™ LT Reusable Underpads
Enhance patient care with high performance barrier

Promote the Ideal Skin Environment


More than 50 percent of older Americans struggle with incontinence.¹ It can increase the risk factors associated with healthcare-acquired pressure injuries, which occur in more than 2.5 million people each year at a cost of \$43,180 per incident.²

1. Moisture

Maintain Comfortable Conditions.

PerforMAX LT creates a dryer, cooler and more comfortable environment for patients by rapidly wicking fluid away from the patient to make the top feel dryer quicker.

How It Works—The Layers of PerforMAX LT

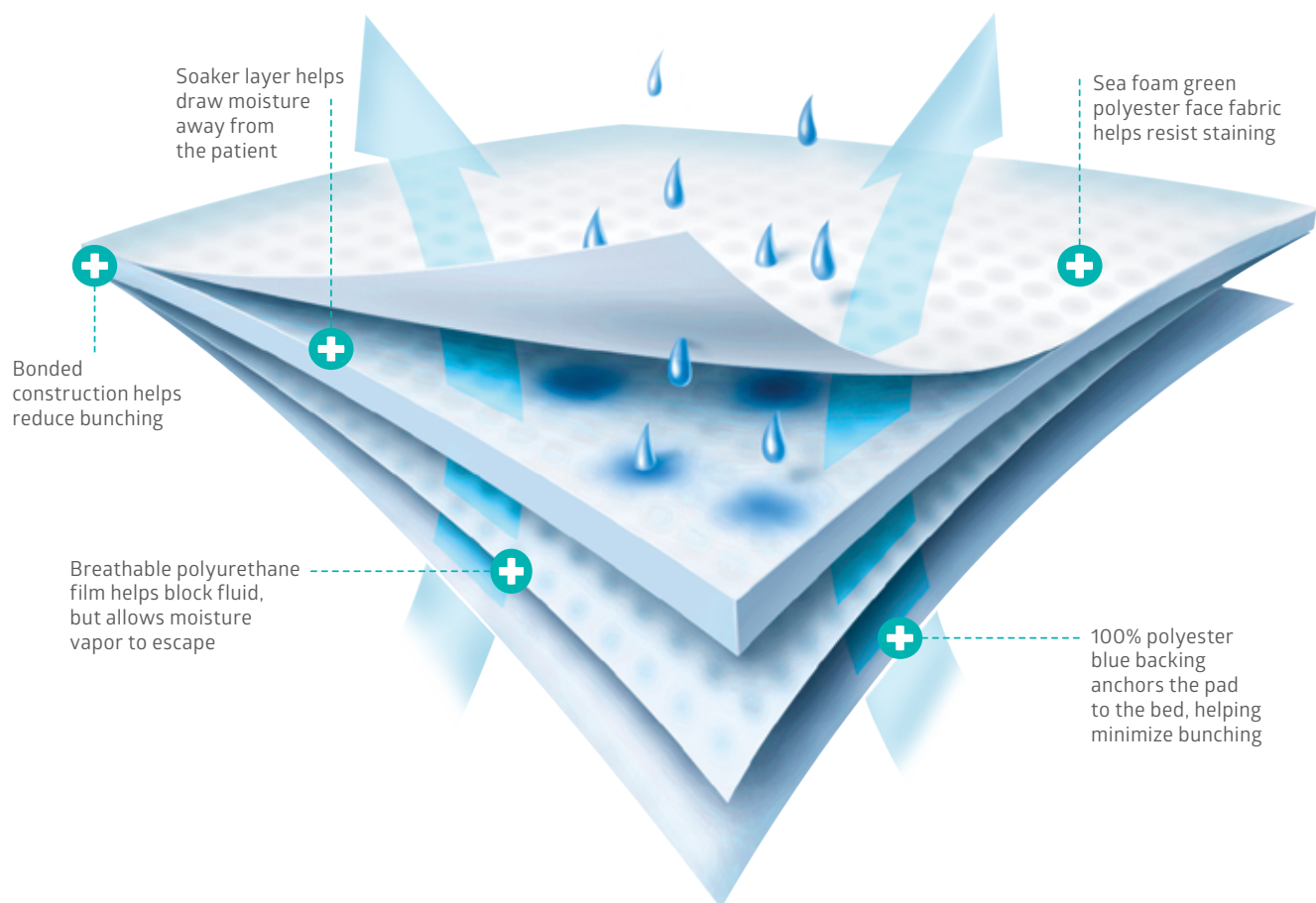

nt. It's easy as 1, 2, 3

But the problem is preventable.


The PerforMAX LT reusable underpad features a cutting-edge moisture-barrier technology to help manage the 3 key factors supporting skin care: moisture, pressure and heat.

2. Pressure


Effective Pressure Management.

PerforMAX LT works with your hospital beds and linens to help support skin care and patient comfort. These pressure maps show pressure distribution is unaffected with the use of the PerforMAX LT pad.⁴ The bonded construction helps eliminate bunching, which can cause added pressure, while the 100% brushed polyester fabric face remains soft against the skin.

Pressure Subject 06—Female 5'3 153 lbs (mmHg – unit of measure for pressure)

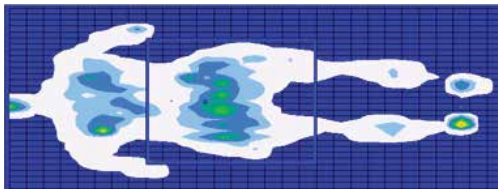

Baseline Average 13.42 mmHg

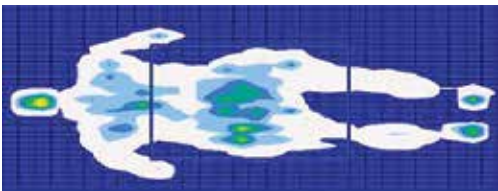

PerforMax Average 11.59 mmHg


3. Heat

Effective Temperature Management.


As these temperature maps show, the use of PerforMAX LT keeps patients nearly as cool as if no pad were used.⁴ Keeping patients cool is crucial to promoting the ideal skin environment.

These temperature maps measure the amount of heat that passes through the pad to a temperature-sensing mat underneath. The higher the temperature, the more heat is moved away from the patient.

Temperature Subject 01—Female 5'6 148 lbs

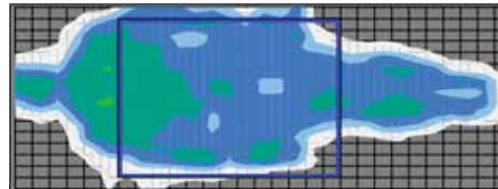

No pad Average 27.21

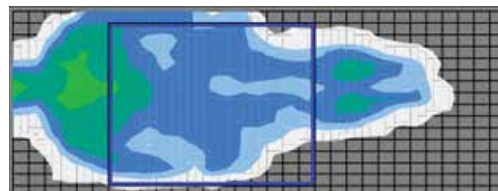

PerforMax Average 26.



More reasons to love PerforMAX LT:

- » Lays flat with minimal bunching or rolling
- » Resists stains so you may get more uses per pad
- » Doesn't interfere with low air loss mattresses³

Ordering Information

Item No.	Description	Pkg
MDTIU7SEFBLU	PerforMAX LT, 100% poly underpad, 34" x 36"	2 dz/cs



Control Odor and Maintain Moisture.

Step up your level of patient care with PerforMAX™ LT SilverTouch®. This specially treated underpad contains all the benefits of PerforMAX LT with the added feature of Silverclear®, a proprietary antimicrobial* solution that reduces odor-causing bacteria.

- » **Prevents Bacterial Growth**—Helps eliminate the source of the odor rather than just masking it
- » **Prolongs Life of Underpad**—Odor control application withstands repeated washings and requires no special laundering

Item No.	Description	Pkg
MDTIU7SOCBLU	PerforMAX LT SilverTouch, 100% poly underpad, 34" x 36"	2 dz/cs







*This product is treated with an antimicrobial agent to protect the product. This product does not protect the user or others from disease-causing bacteria.

To learn more, contact your Medline representative, visit us at Medline.com or call 1-800-396-6996.



Medline Industries, Inc. Three Lakes Drive, Northfield, IL 60093

FOLLOW US    

Medline United States
1-800-MEDLINE (633-5463)
medline.com | info@medline.com

Medline Canada
1-800-396-6996
medline.ca | canada@medline.com

Medline México
01-800-831-0898
medlinemexico.com | mexico@medline.com

References 1. U.S. Centers for Disease Control and Prevention's National Center for Health Statistics. WebMD. <http://www.webmd.com/urinary-incontinence-oab/news/20140625/over-half-of-seniors-plagued-by-incontinence-cdc>. Accessed May 19, 2016. 2. Preventing Pressure Ulcers in Hospitals. AHRQ website. <http://www.ahrq.gov/professionals/systems/hospital/pressureulcertoolkit/putool1.html>. Accessed May 19, 2016. 3. Data on file. 4. Study: Evaluation of Pressure, Temperature and Humidity, Mapping with Reusable Incontinent Pads in Normal Healthy Volunteers, Medline Industries, Inc., June 6, 2012. Data on file.

Some products may not be available for sale in Mexico or Canada. We reserve the right to correct any errors that may occur within this brochure. © 2016 Medline Industries, Inc. All rights reserved. Silverclear is a registered trademark of CTT Group Inc. Corporation Canada. PerforMAX is a trademark, and SilverTouch and Medline are registered trademarks of Medline Industries, Inc. MKT1676693 / LIT374R / 200 / DP / 39