

Xalt reusable surgical gowns

Protection, comfort, durability and sustainability all in one gown



**Perform
with
confidence**

Xalt reusable surgical fabric delivers the ultimate in performance

Our Xalt surgical textiles offer a superb combination of protection, comfort, durability and sustainability. Comprised of a proprietary membrane laminated between two layers of polyester material, Xalt provides your team the level of confidence they need to perform their best.



Durable





Industrial laundry
75 wash cycles




Breathable

Ultimate protection

All products featuring the Xalt fabric and seam seal tape comply with Level 4 regulations throughout their lifespan.

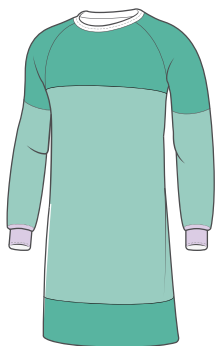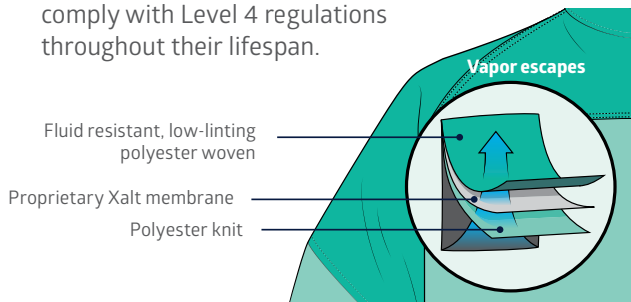

Panel coverage gown

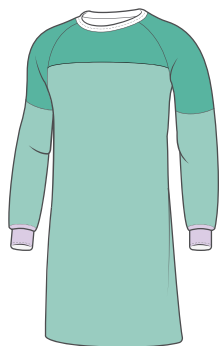

Critical coverage gown



Additional Information

Latex free:	Not made of natural rubber latex
Sterility:	Non-sterile
Traceability:	Lot number on each case of product
Expiration date:	No expiration date
Storage conditions:	General warehouse conditions
Fabric weight:	5.25 OY




Gowns feature a quality control grid that allows you to track usage

Xalt Level 4 surgical gowns


- Highly protective and breathable
- Liquid- and viral-resistant Xalt fabric for critical zones
- Liquid- and static-resistant ResiStat fabric in non-critical zones
- Reusable through 75 wash cycles; low linting fabric



Snaps at collar

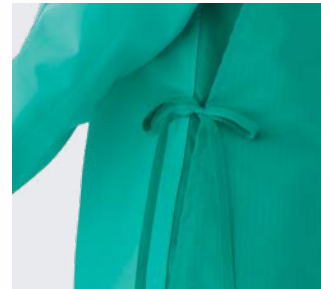

Ties at mid-back

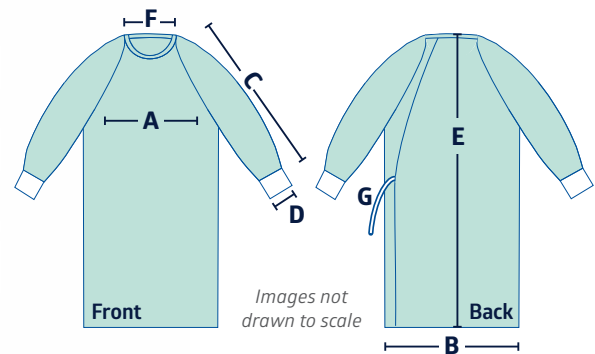

Gown Measurements

	SKU: MDTGXC4JL	MDTGXC4JXL	MDTGXC4JXXL
Size:	Size L	Size XL	Size XXL
Qty:	12 dz/cs	12 dz/cs	12 dz/cs
Gown Color:	Green	Green	Green
Collar Color:	Ceil Blue	White	Green
A Chest	25"	29"	33"
B Sweep	71.125"	80.125"	89.125"
C Sleeve Length	31.125"	34.5"	38.375"
D Sleeve Opening	11.625"	12.25"	12.75"
E CBL	51"	51.875"	52.875"
F Collar	27.125"	29.5"	31.5"
G Cuff	21.5"	21.5"	21.5"

Xalt Level 4 panel coverage gowns




MDTGXP4JXL



Ties at collar



Ties at mid-back

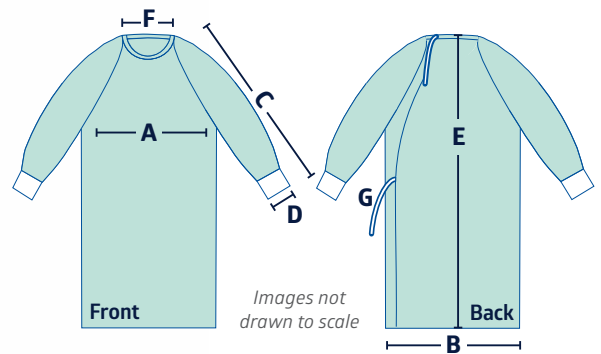

Gown Measurements

	SKU: MDTGXP4JL	MDTGXP4JXL	MDTGXP4JXXL
Size:	Size L	Size XL	Size XXL
Qty:	12 dz/cs	12 dz/cs	12 dz/cs
Gown Color:	Green	Green	Green
Collar Color:	Ceil Blue	White	Green
A Chest	25.25"	27"	29"
B Sweep	63"	67"	72.75"
C Sleeve Length	32"	32.75"	34"
D Sleeve Opening	11"	12"	13"
E CBL	51"	51"	51"
F Collar	27"	29"	32"
G Cuff	21.5"	21.5"	21.5"

Laundering recommendations for Xalt surgical gowns

Normal wash

Operation	Temperature	Wash time
Suds	72°C/162°F	40 min.
Rinse	60°C/140°F	5 min.
Rinse	60°C/140°F	10 min.
Wet Tumble	–	5 min.

DO NOT use softener.
For best results use batch washer.

Bleach

Do not bleach—may cause damage to barrier properties

Normal

Tumble dry: 66°C +/-5°C | 150°F +/-10°F

MAX temperature :71°C/160°F

With final 10 minute cool-down until fabric is dry

Iron

Do not flat-work iron

Prevac Steam Sterilization

Product provided non-sterile

To sterilize product: 134°C/274°F

4 minute exposure

20 minute exhaust

All products:

Must be laundered prior to use

Only wash with like items

Only wash with like colors

AAMI Life Cycle:

Products will pass labeled AAMI Level for up to 75 processings when following above recommendations

DO NOT use fabric softener. Fabric softener inhibits absorbancy/wicking
DO NOT overdry. Overdrying can lead to melting/hardening, sheen, static and permanent creasing

*Surgical materials meet the Standard for Flammability of Clothing Textiles, CPSC 16 CFR Part 1610, Class 1 (formerly CS 191-53).



**Let us know
how we
can help**

Our clinical specialists can analyze your protective apparel usage trends and recommend ways to optimize utilization and purchasing to help you achieve significant savings.

Protect your team today. For more information and samples, contact your Medline Representative or visit [medline.com](https://www.medline.com).



**We make
healthcare
run better™**

Medline Industries, LP
Three Lakes Drive, Northfield, IL 60093
Medline United States | 1-800-MEDLINE (633-5463)
[medline.com](https://www.medline.com) | info@medline.com

Medline Canada
1-800-268-2848 | [medline.ca](https://www.medline.ca) | canada@medline.com
Medline México
01-800-831-0898 | [medlinemexico.com](https://www.medlinemexico.com) | mexico@medline.com

Follow us    